

A bright starburst graphic with multiple rays emanating from a central point, positioned above the 'Rex Point' logo.

Rex Point<sup>®</sup>

STAINLESS STEEL  
SUTURE NEEDLES





## C O N T E N T S

---

04 | TYPE OF NEEDLE

05 | TYPE OF CIRCLE · TYPE OF EYE

06 | 외과용 각침 3/8 CIRCLE(GENERAL CLOSURE) · 외과용 각침 1/2 CIRCLE(GENERAL CLOSURE)

07 | 외과용 각침 STRAIGHT(GENERAL CLOSURE) · 내복용 환침 2 1/2 CIRCLE(ABDOMINAL)

08 | 내장용 환침 1/2 CIRCLE(INTESTINAL) · 내장용 환침 STRAIGHT(INTESTIANL)

09 | 부인과용 환침 1/2 CIRCLE(OB-GYN) · 부인과용 각침 1/2 CIRCLE(VASCULAR)

10 | 안과용 각침 3/8 CIRCLE(OPHTHALMIC) · 혈관용 환침 1/2 CIRCLE(VASCULAR)

11 | 피부용 각침 3/8 CIRCLE(SKIN) · 각막용 각침 3/8 CIRCLE(CORNEAL)

## TYPE OF NEEDLE



### TAPER POINT NEEDLE

THIS TYPE OF NEEDLE PENETRATES, PASSES THROUGH TISSUES BY STRECHING WITHOUT TISSUE CUTTING, NORMALLY THIS NEEDLE MINIMIZES POTENTIAL TEARING OF FASCIA AND IS USED FOR EASILY PENETRATED TISSUES AS LIKE ABDOMINAL VISCERA, PERITONEUM, DURA, SUBCUTANEOUS LAYERS.



### CONVENTIONAL CUTTING NEEDLE

THIS TYPE OF NEEDLE HAS 3 CUTTING EDGES AND DESIGNED FOR TISSUE THAT IS TOUGH TO PENETRATE AS LIKE SKIN, TENDON SHEATHS, ORAL MUCOSA.



### REVERSE CUTTING NEEDLE

SAME WITH ABOVE CONVENTIONAL TYPE AND THESE NEEDLES ARE ALSO HELPFUL IN COSMETIC AND OPHTHALMIC SURGERY, CAUSING IT'S MINIMAL TRAUMA.




### SPATULA CUTTING NEEDLE


THIS TYPE OF NEEDLE DESIGNED TO REDUCED TISSUE INJURY AND ALLOW MAXIMUM EASE OF PENETRATION AND CONTROL AS IT PASS BETWEEN AND THROUGH TISSUE LAYERS.

## TYPE OF CIRCLE

135° (3/8)

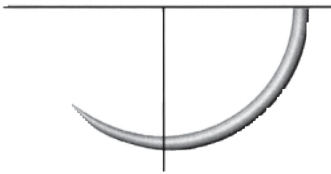

180° (1/2)

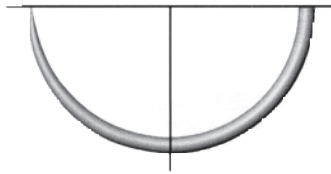


STRAIGHT



## TYPE OF EYE

REGURAL EYE

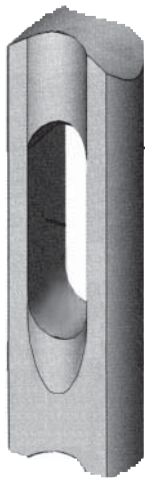

SPRING EYE






외과용 각침 3/8 CIRCLE  
(GENERAL CLOSURE)

102 : REGULAR EYE


101 : SPRING EYE

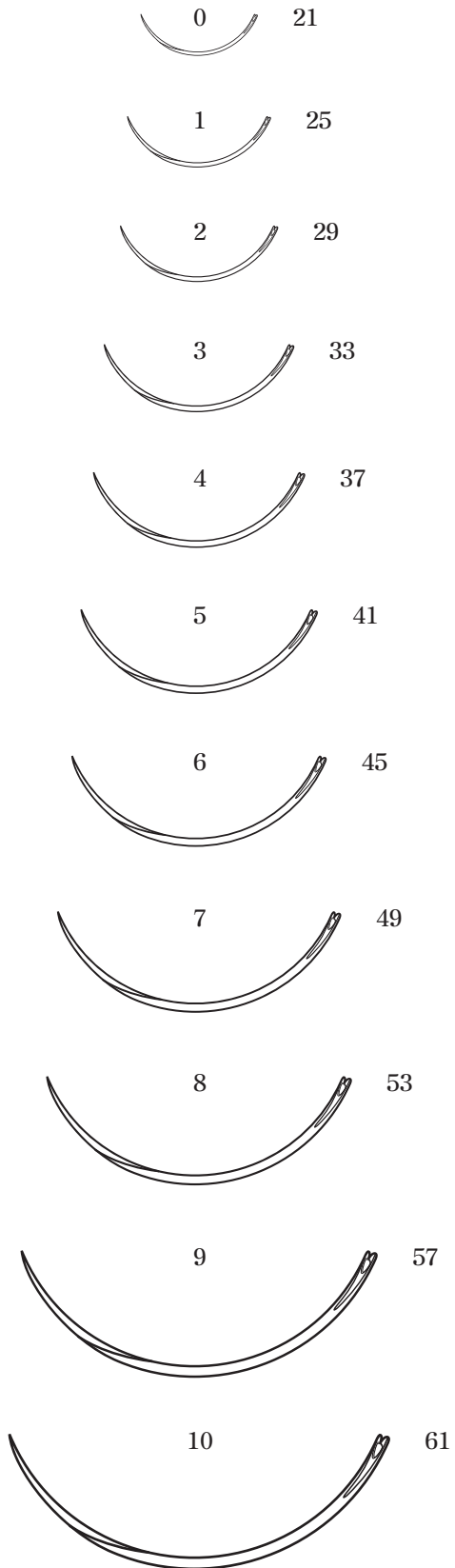

외과용 각침 1/2 CIRCLE  
(GENERAL CLOSURE)

106 : REGULAR EYE

105 : SPRING EYE

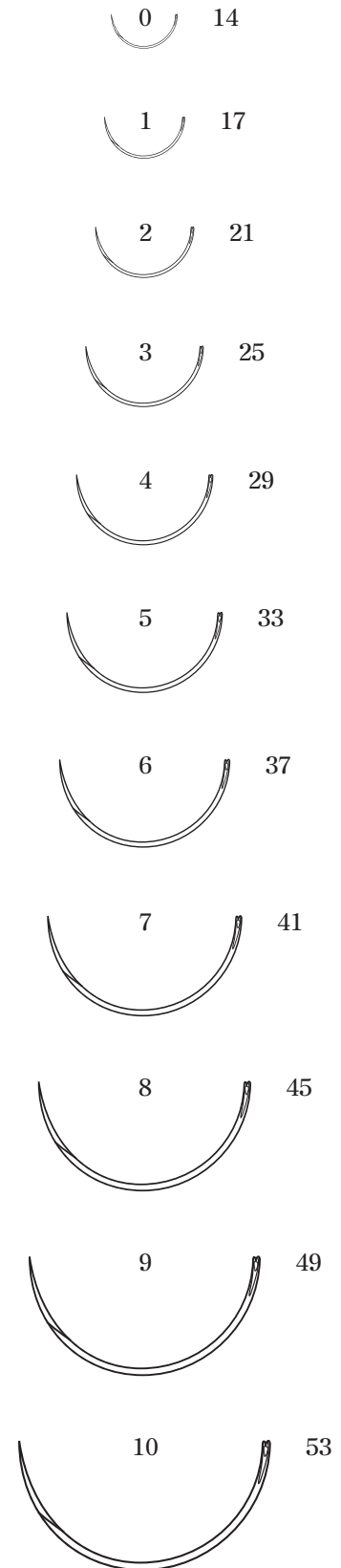




외과용 각침 STRAIGHT  
(GENERAL CLOSURE)

162 : REGULAR EYE


161 : SPRING EYE

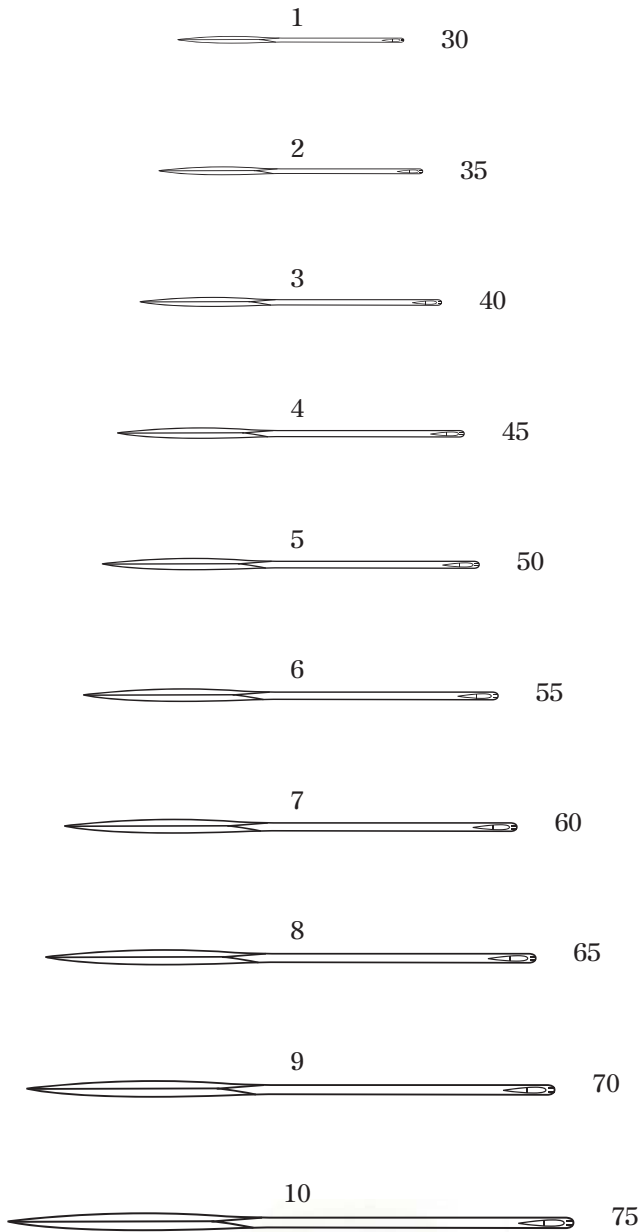

내복용 환침 1/2 CIRCLE  
(ABDOMINAL)

114 : REGULAR EYE

105R : SPRING EYE






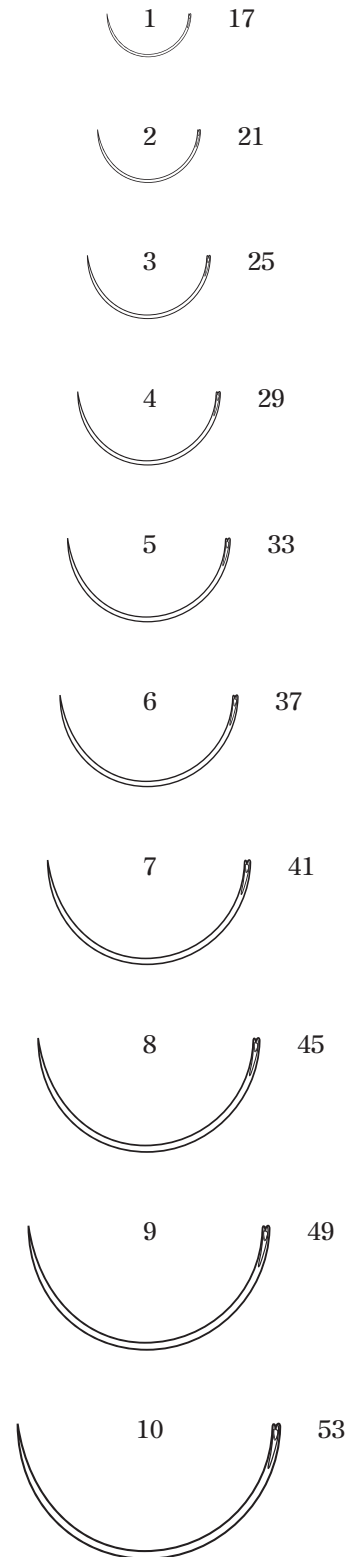




내장용 환침 1/2 CIRCLE  
(INTESTINAL)

136 : REGULAR EYE








135 : SPRING EYE

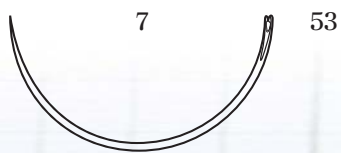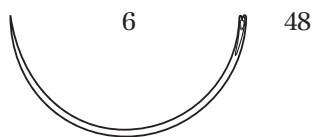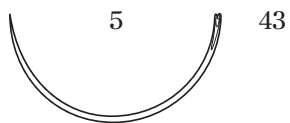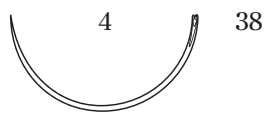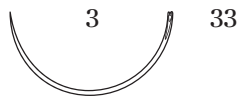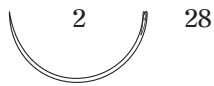

내장용 환침 STRAIGHT  
(INTESTINAL)

134 : REGULAR EYE

133 : SPRING EYE

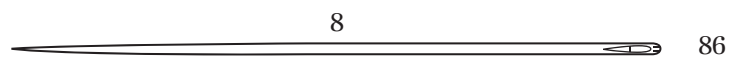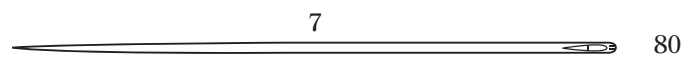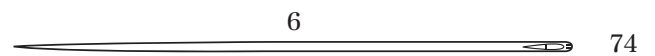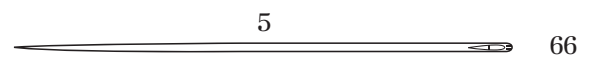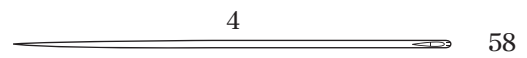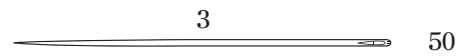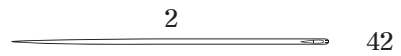











부인과용 환침 1/2 CIRCLE  
(OB-GYN)

144 : REGULAR EYE







143 : SPRING EYE

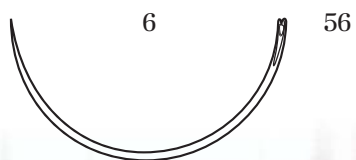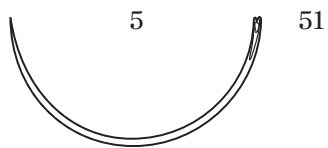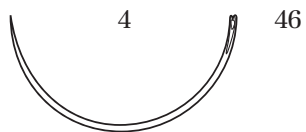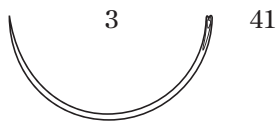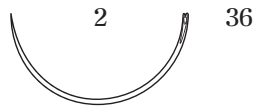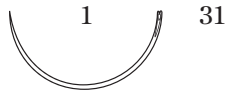

부인과용 각침 1/2 CIRCLE  
(OB-GYN)

142 : REGULAR EYE

141 : SPRING EYE

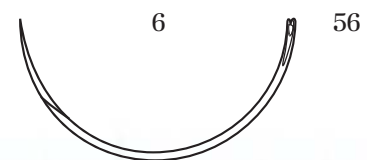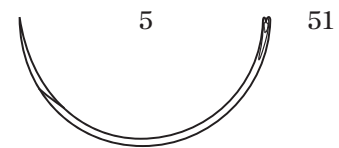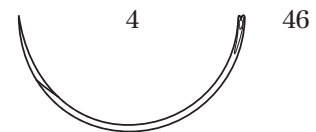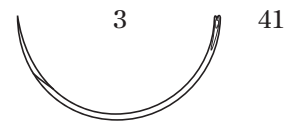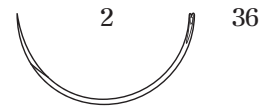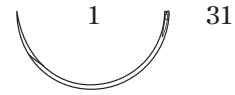



안과용 각침 3/8 CIRCLE  
(OPHTHALMIC)

124 : REGULAR EYE

123 : SPRING EYE

000 9

00 11

0 13

1 15

2 18

3 21



혈관용 환침 1/2 CIRCLE  
(VASCULAR)

24 : REGULAR EYE

22 : SPRING EYE

0 11

1 14

2 19

3 24

4 29







피부용 각침 3/8 CIRCLE  
(SKIN)

104 : REGULAR EYE

103 : SPRING EYE

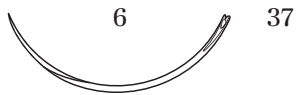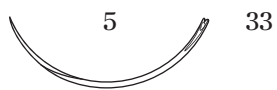

각막용 각침 3/8 CIRCLE  
(CORNEAL)

128 : REGULAR EYE

127 : SPRING EYE





356-3, SAMNAK-DONG, SASANG-GU, BUSAN, KOREA TEL : 81-51-305-8016 FAX : 81-51-305-8021  
E-mail : sales@sm-eng.net www.sm-eng.net